

#whdentalwerk



video.wh.com

Now at  
your dealer  
or [wh.com](http://wh.com)



# The formula for a long instrument life



The more  
intelligent formula

# Prolong the lifespan of your instruments

## For every application

W&H Service Oil F1 is a product of W&H research. It has been developed hand in hand with modern W&H instrument technology. This oil gives you the confidence of knowing that you are providing your W&H instruments with the best possible care:

- › Turbines
- › Straight and contra-angle handpieces
- › Surgical instruments
- › Air scalers
- › Air motors



W&H dental drives achieve top speeds of up to 390,000 r.p.m. F1 provides excellent lubrication and internal cleaning thereby ensuring fault-free operation of the instruments.



**Perfect for manual maintenance**  
W&H Service Oil F1 in the MD-400 spray can.

# With the test champion oil for all standard drives


## The more intelligent Service Oil

W&H Service Oil F1 is an all synthetic, special oil of maximum purity, which has been specially developed for the demanding and extreme requirements of dental instruments. Compared with other well-known brand oils, the Service Oil F1 was subjected to a series of rigorous tests.

The Result: Test winner in all significant performance features. The test result: The better the oil, the lower the wear. (Indentation diameter of friction WK [mm] determined from TU-Wien according to the SRV-test (translatory oscillation-test) standard DIN 51 834.)

## Short lubrication intervals – long service life

Intensive maintenance will considerably increase the life of your instruments.

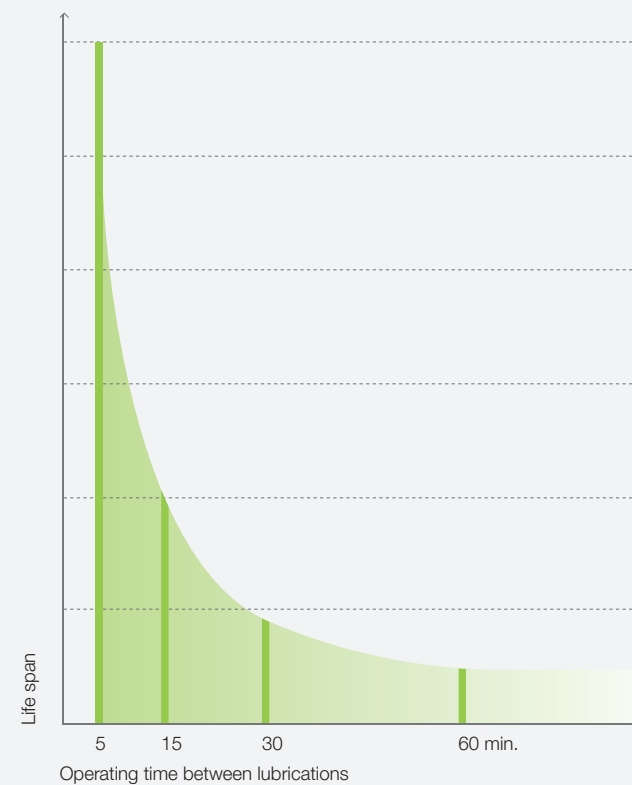

## Sterilization resistant 135°C

F1 is of course fully sterilizable, absolutely resistant to hydrolysis, stable against oxidation and is shown to be highly compatible with disinfectants.

## Physiologically safe

F1 contains no toxic constituents, is physiologically harmless and environmentally compatible. The substances it contains comply with the applicable international regulations.

## High material compatibility

F1 also safeguards the durability of plastic parts, coatings, seals, glass etc.

## Enormously economical

Your instruments will thank you for regular maintenance with F1 by requiring fewer repairs. The advantages: considerable savings in cost and increased profitability of your W&H instruments.

# Intelligent oil + short maintenance intervals = longer lifespan

**Oil spray can** .....  
Type MD-400  
Contents 400 ml  
also available in  
packs of 6




Order your Service Oil F1 through your dental supplier.  
Suitable spray caps can be found on the Internet at **wh.com**

## Distribution:

**W&H Dentalwerk Bürmoos GmbH**  
Ignaz-Glaser-Straße 53, Postfach 1  
51111 Bürmoos, **Austria**  
t +43 6274 6236-0  
f +43 6274 6236-55  
office@wh.com  
**wh.com**

**W&H (UK) Limited**  
6 Stroud Wood Business Centre  
Park Street, St Albans, Hertfordshire AL2 2NJ, **United Kingdom**  
t +44 1727 874990  
f +44 1727 874628  
office.uk@wh.com  
**wh.com**

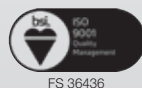

Pictures are for illustrative purposes only. Additional equipment and accessories shown are not included as standard.

20337 UEN Rev. 004 / 27.02.2024  
Subject to alterations