

#whdentalwerk



video.wh.com

Now at
your dealer
or **wh.com**



Reliable. Efficient. Precise.

Turbines



alegra

TRIPLE SPRAY

for optimal cooling of the treatment site from all sides.

OPTIMAL LED POSITION

for optimal view and illumination of the treatment site.

HIGH TORQUE

for fast removal of tooth structure.

PRECISION MANUFACTURE

for vibration-free and quiet operation.

SMALL TURBINE HEAD

for optimal view of the treatment site.

PUSH-BUTTON CHUCKING SYSTEM

for simple, rapid bur changes with low actuating force.

CERAMIC BALL BEARINGS

for extremely quiet operation and a long service life.

WATER FILTER

for perfect spray function.

INNOVATIVE GENERATOR SYSTEM

in the Roto Quick quick coupling for use with Alegra light turbines.

Long service life and extremely quiet operation


Discover the new Alegra turbines with extremely quiet operation for daily use in your practice. Designed to meet your needs with state-of-the-art technology. Made in Austria: For your best treatment results.

Made in AUSTRIA



Long service life

Ceramic ball bearings guarantee minimum wear. Particularly with optimum care and maintenance!

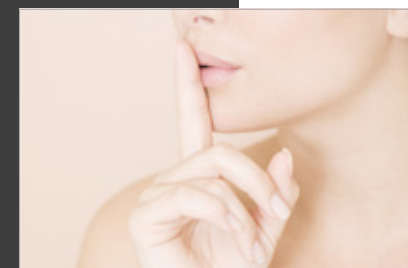


Quiet

Alegra turbines offer top quality, vibration-free, quiet operation thanks to precision manufacture.



Replaceable water filter

Prevents dirt particles from entering the turbine.

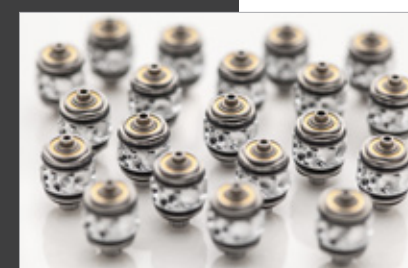

Your reliable partner

This combination makes W&H turbines a reliable partner when working on your patients.

Up to the minute


W&H is a dynamic and rapidly expanding company. We guarantee reliable support and outstanding service. Our recipe for success: innovative technologies, high quality material concepts, stringent quality criteria and many years of experience in the development of precision dental instruments.





Maximum efficiency and safety during treatment

The Alegria range of turbines offers high-end technology for top preparation results! Experience maximum efficiency during use combined with maximum safety: for you and your patients.

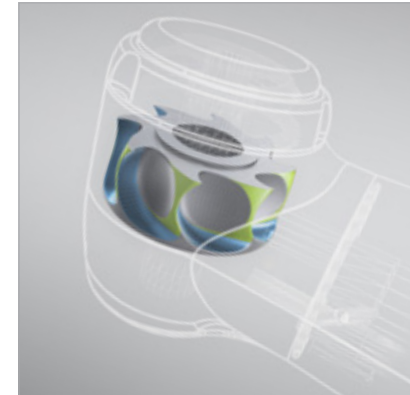

Full steam ahead with the Power Blade rotor

The optimized blade geometry allows higher rotor efficiency. For maximum torque and optimal power throughout preparation.



Precision cooling

The triple spray during treatment cools the treatment site from all sides and guarantees safety during treatment.



Tailor-made turbines

From hard-to-reach areas to fast removal of tooth structure: The right head size for every application.

Hygiene head system

The circulating air in the turbine head ensures no aerosol particles are sucked into the turbine as the rotor slows down.



Simple bur changing


Optimum balance between easy operation of the push-button and maximum retention force of the chucking system. No tools required!

Hygiene & maintenance of your turbines

“The better the care, the longer the pleasure.” Regular cleaning and care of your instruments protects them against premature wear and unnecessary service costs.



Short lubrication intervals – long service life

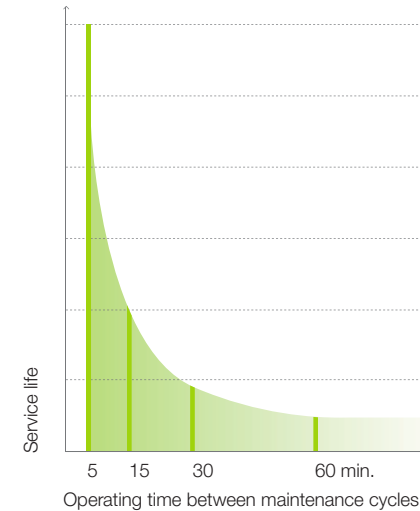

For every application

W&H Service Oil F1 is a product of W&H research. It was developed hand in hand with the leading edge W&H drive technology. This oil gives you the security of knowing that you are providing your W&H products with the best possible lubrication.

Perfect for manual lubrication

W&H Service Oil F1 in the spray can cleans and lubricates your turbines to guarantee troublefree use.

Maintenance at the press of a button

The range of W&H maintenance units offers you optimal lubrication for your turbines. And all in record time!

Optimal view and independent LED

The integration of a generator in the Roto Quick quick coupling gives you a low-cost option of adding an LED turbine to your selection of Alegria turbines. The advantages for you: Flexibility in use and in hygienic preparation.



Excellent ergonomics

Fatigue-free working thanks to easy 360° rotation of the turbine on the coupling, low overall weight and slender design.

Further advantages at a glance

- › Integrated resuction stop: Prevents resuction of aerosol particles into the spray channel
- › Component replacement on site by the user
- › Thermo washer disinfectable
- › Sterilizable up to 135 °C
- › 12 months warranty



LED turbine + quick coupling with generator = independent light

A generator in the Roto Quick quick coupling enables the use of light turbines on the turbine hose without a power supply.



Optimal view

The position of the LED in the turbine head guarantees optimal view and illumination of the treatment site.

Technical data



Turbines with LED: Couplings with generator: Turbines without light:	TE-97 LQ* RQ-53 / RQ-54 TE-97 / TE-97 RM / TE-97 BC	TE-98 LQ* RQ-53 / RQ-54 TE-98 / TE-98 RM / TE-98 BC	TE-95 RM / TE-95 BC
Light (LED version):	LED and generator		–
Head size:	Ø 10.4 mm	Ø 12.2 mm	
Power**:	16 W	18 W	
Speed (min ⁻¹):	390,000 ± 30,000	330,000 ± 30,000	
Rotary instruments:	FG bur, in acc. with ISO 1797-1, Ø 1.6 mm		
Max. working part diameter:	2 mm		
Max. length:	21 mm	25 mm	
Weight:			
› Turbines with LED	46 g	48 g	
› Couplings with generator	26 g / 20 g	26 g / 20 g	
› Turbines without light	46 g / 47 g / 47 g	48 g / 50 g / 50 g	49 g / 49 g
Operating pressure:	2.2–2.8 bar (recommended: 2.5 bar) (32–40 PSI)		
Ball bearings:	Ceramic		Steel
Spray:	Triple		Single

* TE-97 LQ and TE-98 LQ turbines with LED are intended for use with W&H RQ-53 and RQ-54 Roto Quick couplings only.

** Performance value applies at a rated pressure of 2.5 bar for turbines without light, of 2.8 bar for turbines with light and exhaust air pressure of max. 0.2 bar. When operating with the maximum permissible operating pressure of 2.8 bar, it is possible to reach 16 watt for TE-97, TE-97 RM, TE-97 BC and 18 watt for TE-98, TE-98 RM, TE-98 BC.

Attention: Not all products mentioned in this brochure are available in all countries. Please contact your local W&H representative for details.

Manufacturer:

W&H Dentalwerk Bürmoos GmbH

Ignaz-Glaser-Straße 53, Postfach 1

5111 Bürmoos, **Austria**

t +43 6274 6236-0

f +43 6274 6236-55

office@wh.com

wh.com

Pictures are for illustrative purposes only. Additional equipment and accessories shown are not included as standard.

20415 AEN Rev. 008 / 03.07.2026
Subject to alterations