

#niwop  
video.wh.com

Now at  
your dealer  
or [wh.com](http://wh.com)



# The new generation of rotary polishing



**proxeo** TWIST **LatchShort**  
Polishing System

## Tiny but mighty

The Proxeo Twist LatchShort Polishing System is a particularly smart solution for Prophy polishing: The small head with the optimally matched Prophy Cups significantly reduces the working height. You gain as much as four millimetres compared to standard latch systems.



### Access all areas


The molar region is more easily accessed from a distal to a mesial direction.

### Comfortable and simple

The easier the access the more comfortable the treatment.

### Particularly suitable for narrow areas

The molars in the tuber region are easier to reach.

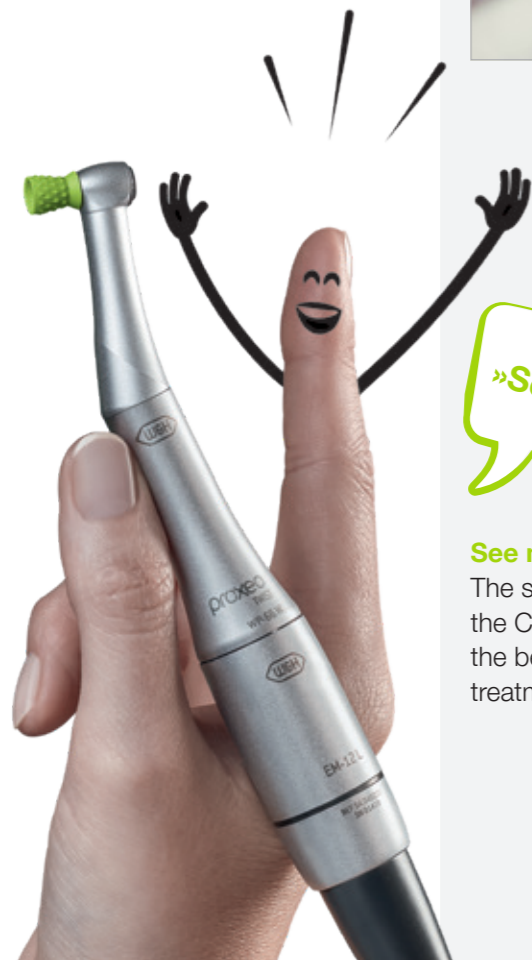

»Sooo much space«

### See more


The smaller the head, the shorter the Cup and the thinner the neck, the better the view of the treatment site.

## Perfect harmony

Smaller. Finer. Innovative: The combination of the slender W&H Proxeo Twist LatchShort contra-angle handpiece with its small head and the Prophy Cups with a shorter shaft make access easy, particularly in narrow areas and to the posterior molars.




Green contra-angle handpiece vs LatchShort system

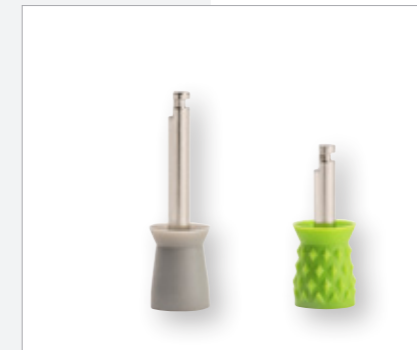

### Measurably more space

The advantage of the new Proxeo Twist LatchShort system compared to the green contra-angle handpiece: More space for polishing with the shorter shaft.

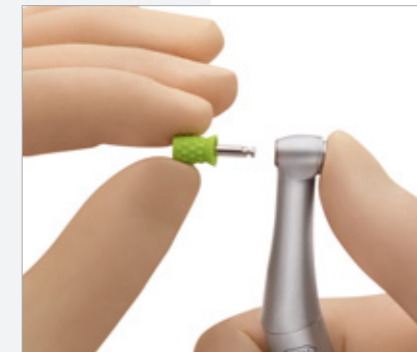

### Change at the push of a button

With the practical push-button chuck system, you can change the Prophy Cups in seconds.



### Perfect sealing

Saliva and polishing paste can no longer infiltrate the head so easily. This supports the long service life of the contra-angle handpiece.



## Inspired by Prophy professionals

In addition to the Proxeo Twist contra-angle handpiece, W&H is now also offering Prophy Cups and Brushes. They were designed in consultation with Prophy users. Great emphasis was placed on features such as optimal adaptation, simple paste pickup and distribution as well as gentle cleaning up into the sulcus.

## One-stop shop

The new Proxeo Twist LatchShort Prophy Cups in various hardness levels and Brushes are now available for optimal removal of biofilm and discolouration. All Proxeo Twist LatchShort Prophy Cups have a shorter shaft for an even lower working height.




**Smooth surface at the end of the Cup** for gentle polishing up into the sulcus.



**Special nub shape** designed to prevent the paste from spraying.



**Special internal fins** for easy pickup and targeted application of paste.

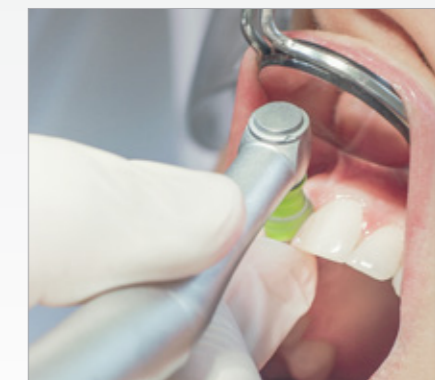

**Nubs on the outside** polish the interdental space at the same time.

## Irresistible advantages for rotary polishing

### Slimmer. Finer. Innovative.

The new Proxeo Twist LatchShort Polishing System with a working height that is up to 4 mm shorter thanks to the Proxeo Twist Prophy Cups and Brushes: for better access, more space and an optimal view of the treatment area.



## Maximum protection and comfort



### Slim and flexible

The contra-angle handpiece easily fits into the hand with its ergonomic design and slim neck. The combination with the W&H EM-12 L electric motor adds to its ergonomic properties.



### Hygiene

The smooth exterior design makes the sterilizable and thermo washer disinfectable Proxeo Twist LatchShort contra-angle handpiece particularly easy to clean. The Prophy Cups and Brushes are disposable products and are not intended for reuse.

### Technical data

Type	WP-66 W	CE 0297	W&H	ISO 1350
Max. speed:	10,000 rpm			
Transmission ratio:	4:1			
Recommended speed for Prophy Cups and Prophy Brushes	1,250 – 2,500 rpm			
Attachments that can be used:	W&H Proxeo Twist LatchShort Prophy Cups and Prophy Brushes			

### Prophy Cups and Brushes



Not for reuse



Latex-free



## Melanie Ruckenthaler

W&H product manager


»The best of both worlds: the unique combination of a low working height and convenient push-button chuck system combines the advantages of screw and latch systems: easier access, better view and maximum comfort on the one hand and the ease of changing polishing Cups and Brushes on the other.«



## Monika Riedl

dental hygienist, Munich

»The Proxeo Twist LatchShort is the first system with a low working height and a push-button. With the working height as much as 4 mm lower than standard latch systems, there is more space and a perfect view of the treatment site.«

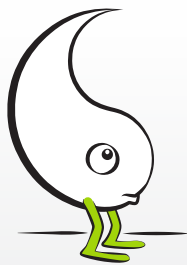

Implantology  
without **Periodontology** is  
like **Yin** without **Yang**.\*



**#NIWOP:** The systematic and evidence-based workflow

### W&H Dentalwerk Bürmoos GmbH

Ignaz-Glaser-Straße 53, Postfach 1

5111 Bürmoos, **Austria**

**t** +43 6274 6236-0

**f** +43 6274 6236-55

office@wh.com

**wh.com**