

#whdentalwerk



video.wh.com

Now at
your dealer
or wh.com



The new solution in **the Zygoma** **application**



Unique, easy and safe

The specially for the insertion of Zygoma implants developed 20:1 handpiece facilitates the work of the surgeon. The angulation is unique and allows a much better view of the treatment site. So the surgical procedure is also possible to difficult accessible positions in the zygomatic bone.

Easy access to treatment site

The angle of the handpiece allows optimal access for difficult interventions in the area of the zygomatic bone.



High application security

The world's only 20:1 handpiece with contra-angle chucking system ensures application security at high concentricity.



Fatigue-free working

The ergonomic design of the handpiece ensures long and fatigue-free working.




Efficient cooling

The spray exits near the instrument's area of activity and cools the operating site extremely efficiently.



Technical data Zygoma

Type	SZ-75
Transmission ratio:	20:1
Coupling system:	ISO 3964 (DIN 13940)
External cooling system:	single
Chucking system:	twist chuck system
Rotary instruments:	Zygoma burs, surgical burs and cutters with contra-angle shank Ø 2,35 mm
Maximum drive speed:	50,000 min ⁻¹
Can be dismantled:	yes
Indications, e.g.:	Implantation of Zygoma implants osteotomy

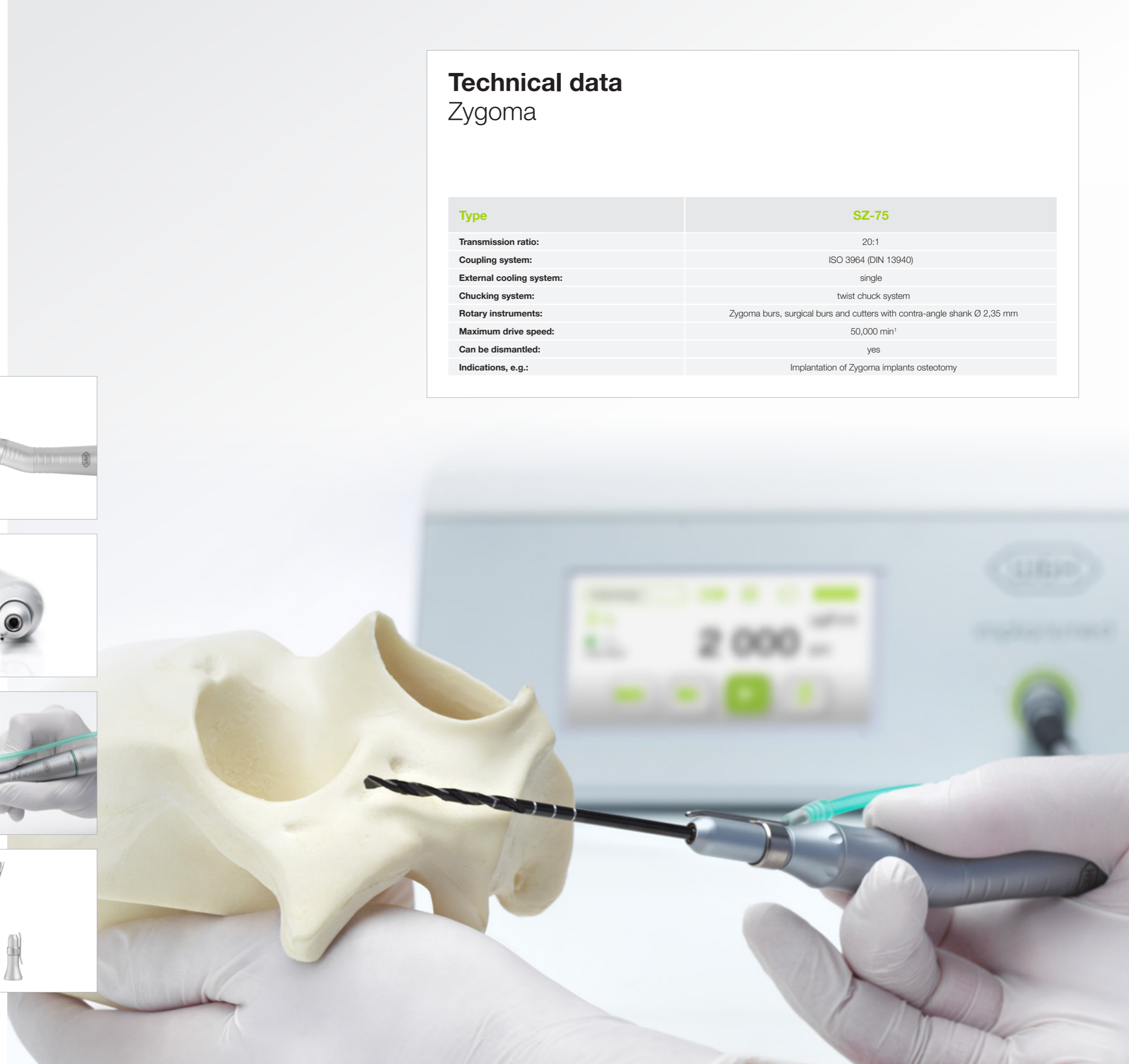



Dr Armando Lopes

Specialist in Oral Surgery & Implantology
MALO CLINIC, Portugal

“W&H offers the best handpiece for the insertion of Zygoma implants: The unique combination of the specially angled handpiece, the efficient spray cooling and a secure chucking system allows maximum comfort. For doctor and patient.”



Dr Ana Ferro

Specialist in Periodontology
MALO CLINIC, Portugal

“Beautifully machined for quality and balance. W&H sets new standards in the Zygoma application. Incredibly easy to use with Implantmed and fatigue-free working thanks to the ergonomic design of the handpiece are guaranteed.”



W&H Dentalwerk Bürmoos GmbH

Ignaz-Glaser-Straße 53, Postfach 1

5111 Bürmoos, **Austria**

t +43 6274 6236-0

f +43 6274 6236-55

office@wh.com

wh.com